

# CATALOG

## CHOUKROUN PRF™ SYSTEM



**PROCESS FOR PRF™**  
CHOUKROUN

49, rue Gioffredo • 06000 NICE • FRANCE  
www.a-prf.com • info@a-prf.com  
00.33.493.855.890

# THE EQUIPMENT

## ▶ DUO Quattro Centrifuge:

Customizable centrifuge !




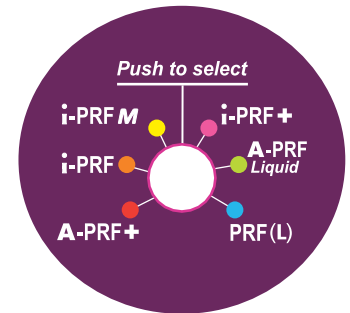
❄ Winter 🍂 Autumn 🌸 Spring ☀ Summer

Change the panels according to the seasons !

The new centrifuge DUO Quattro has a «push» button that affords you to select between 6 different PRF protocols. The values of these settings are preset but not locked (you can set other values if needed) .

**Additional position** : There is a 7<sup>th</sup> position after PRF (L) program, it's a manual settings position. You can set and memorise the values you want.

Simply push the button to select the desired mode of operation, the corresponding LED will light on. The machine is ready. Simply press the button  to start the centrifugation.



**6 PROTOCOLS**

## ▶ PRF BoX :



« PRF Box » allows to get membranes always hydrated and of constant thickness, but also to recover the exsudate rich in proteins: Vitronectin and Fibronectin. You can also produce “plugs” of PRF.

**NEW** : PRF box with a line hole for the punching of the membrane.





## The instruments :



### PRF Scissors

Allows to separate easily the fibrin clot from the red cells and to cut membranes into small pieces.



### PRF Forceps

Is used to remove and manipulate membranes.



### PRF Forceps giraffe

The forceps giraffe is intended to tunneling the PRF membranes for surgical soft tissues procedures: introduce the forceps in the tunnel, open the forceps, block the PRF membranes and pull to place the membranes in the correct position.



### PRF Mini Tray

- Used to prepare membranes and too:
- Separate the fibrin clot from the red blood cells
  - Cut the membranes
  - Fold the membranes
  - Make «punches»



### PRF Bowl

Is used to mix the membranes cut into small pieces with biomaterials. It can also be used for the cutting: Put the membrane into the bowl and cut it with PRF scissors.



# The instruments :

1  
In the Sinus: this plate is used to retrieve the PRF membrane (always) folded twice and slide it into the sinus to the desired location easily. Introduce the pointed end of the pad 1 in the sinus window and drag the PRF membrane on the Schneiderian membrane. Same process for a collagen membrane.

## Pad PRF

2  
Is used to position the PRF membranes easily, anywhere, using the pad as a cake server. Drag the PRF membrane on Pad 2 and place it on the desired surgical area. Same process for a collagen membrane...



## Compactor Small

Allows to compact and tamp the bone filling materials in the recipient sites (alveoli, sinus, goshawks implants).



## Compactor Big

Allows to compact and tamp the bone filling materials in the recipient sites (alveoli, sinus, goshawks implants).



## Double Spoon

Allows to collect and handle the bone filling materials in the recipient sites (alveoli, sinus, goshawks implants).

▶ The instruments :

## Tube-holder

Is used to hold the tubes and let them rest after the blood drawing



## Tourniquet

It must be tightened to the patient's forearm 10 cm above the site of venipuncture.



## Sticker

Is used to maintain in place the collector during the blood drawing



▶ Tubes and Blood Collectors :



Clot Membrane

**A-PRF™ +**



Liquid Injectable

**i-PRF™**



**Blood Collector**

# POLYSTERIBOX®



The polysteribox® is a reusable container for sterilization, transport and storage of your instruments. This box is made of «DURADEX» (PPSU), This material allows higher stability, durability and resistance. (up to 150°C)

Filter duration **5 years** or **1250 sterilizations**  
Duration of sterilization: **6 months**

## Cleaning:

Clean with an alkaline detergent and rinse with water.

## Sterilization:

- ✓ Autoclave sterilization up to 134°C.
- ✓ Chemical disinfection at a maximum temperature of 60°C
- ✓ Thermal disinfection at a maximum temperature of 93°C

**⚠ Do not use drying aids! Risk of damage! The material becomes opaque.**

The polysteribox® L can contain all the PRF kit for an optimal and practical sterilization.



**Caution:** Remove the crusher from the PRF BoX™ before sterilization




Arrange all Instruments in the Polysteribox®



Close, place the sterilization label. Then proceed with the sterilization. After that seal with the safety seals

## Delivered with :

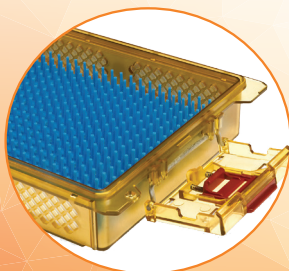

Silicone positioning mat



Sterilization seals



Safety seals

## Optional :



Strainer



49, rue Gioffredo • 06000 NICE • FRANCE  
[www.a-prf.com](http://www.a-prf.com) • [info@a-prf.com](mailto:info@a-prf.com)  
00.33.493.855.890